
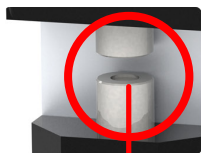


# PORTABLE DIGITAL TABLET HARDNESS TESTER

298/DGP

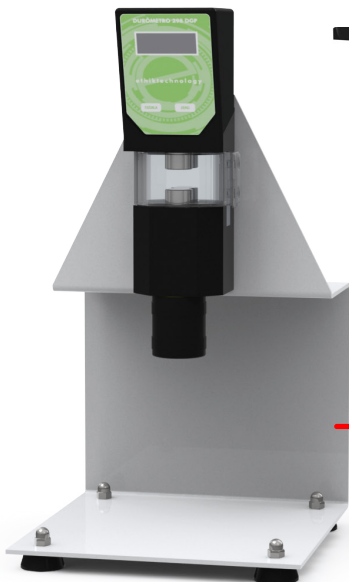

jaw depressed base for  
hardness tests of pel-  
lets, animal food, etc.

Fixation support for  
workbench use.

## Product description

- Equipment developed to perform hardness measuring on pills, fertilizers, candies, animal food, pellets, etc
- Portable and easy handling.

## Technical Characteristics

- Digital microprocessor with liquid crystal viewer (LCD)
- Applied force visualization in real time until tablet's breaking point
- Aluminum body and knob, polycarbonate structure and translucent acrylic sample protector
- Clear test area for easy cleaning
- Automatic shut down after 2 minutes of stand by.

## OPTIONAL



- Fastening rack, depressed area for hardness tests in fertilizers, pellets, animal food and others
- This rack allows the equipment usage on workbench and depressed jaw base assures precision on measurement of animal food, pellets and fertilizers.

## INCLUDED ACCESSORIES AND COMPONENTS



- Operation manual
- Calibration Certificate – traceable RBC
- Four (04) 1,5V batteries

## Technical Specifications

Physical tablet characteristics	Maximum diameter 25mm
Hardness resolution	0,1
Precision	0,1
Scale	4 to 20mA
Hardness measurement unit (manual select)	Kgf; N or Lb
Protection grade	IP 20
Electrical input	4 x 1,5V batteries